


Bottle Cap with Integral Probe Clamp and Cannular for Validation Thermocouples

Product Code: GL45-1





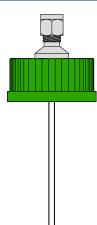

Description:

When using liquid containers fitted with a validation thermocouple in a steriliser, it is often difficult to ensure that the thermocouple hot junction is maintained centrally and at the desired depth within the container during the full sterilising cycle.

The GL45-1 makes this a simple operation by having a permanently fitted integral cannular. The cannular is made to the required length and allows the thermocouple hot junction to be held just below the cannular opening.

A clamp is permanently fitted to the cap which can be supplied as standard or with a vent hole if required.

Specification:	
Moulded polypropylene	
DIN thread GL45	
-40 to +140°C steam.	
Round glass Duran® bottles of 100ml, 250ml, 500ml and 1000ml sizes.	
Square Duran® bottles of 250ml, 500ml and 1000ml.	
Square polycarbonate* Corning® bottles of 150ml, 250ml, 500ml and 1000ml.	
316/316L stainless steel	
In Teflon®.	
In Teflon® tube, 6.0mm OD x 3.8mm ID, length cut to requirement. Stainless steel cannulars also available.	

Duran® is a registered trade mark of Schott Duran Produktions GmbH & Co. KG. Corning® is a registered trade mark of Corning Inc. Teflon® is a registered trade mark of the Dupont Corporation.

Datasheet: 027

Thermal Detection Limited Unit 6, Primrose Hill Industrial Estate Orde Wingate Way Stockton-on-Tees, England, TS19 0GA



Fax: +44 (0) 1642 618307
Tel: +44 (0) 1642 602878
Email: sales@thermal-detection.com
Web: www.thermal-detection.com

© Thermal Detection Ltd 2016



Bottle Cap with Integral Probe Clamp and Cannular for Validation Thermocouples

Validation Accessories:



Pre-cut Validation Thermocouples



Validation Thermocouple Looms



GL45-1 Bottle Cap Gland for Thermocouple



Thermocouple Entry Gland TCEG




Y Piece Extension



Tee Piece Extension



Clamps, Tri-clover Blanks, Gaskets

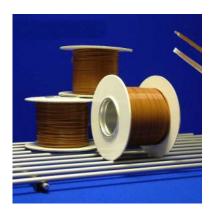

Type T Thermocouple Wire



Ampoule Holders

Datasheet: 027

Thermal Detection Limited Unit 6, Primrose Hill Industrial Estate Orde Wingate Way Stockton-on-Tees, England, TS19 0GA



Fax: +44 (0) 1642 618307 Tel: +44 (0) 1642 602878 Email: sales@thermal-detection.com

Web: <u>www.thermal-detection.com</u>
© Thermal Detection Ltd 2016